EDAQ

e-corder 821 (Model ED821)



- Turns your PC or Macintosh into a precision data recorder
- No programming required plug and play installation
- Powerful on-line and off-line analysis
- Continuous recording up to 200 000 samples/s
- Independently selectable input gains ± 2 mV to ± 10 V
- 16 bit A/D resolution
- Bipolar output for pulse and waveform generation
- Trigger input TTL or contact closure
- USB 2.0 and 1.1 compliant

Description

The e-corder 821 is a high performance eight channel data recording and analysis system. It records analog data from a wide variety of transducers and instruments, along with eDAQ Amps and Pods. It features eight programmable input amplifiers, four of which can be configured for differential input. The system can record 16 bit resolution data directly to your PC at speeds of up to 200 000 samples/s via a USB 2.0 compliant interface. A built-in software controlled analog output, provides basic pulse and waveform generation.

Computer System Compatibility

e-corder is compatible with the following computer environments.

Windows 2000 or XP or later. Systems with a minimum of 128 MB RAM and USB port.

MacOS X or later. Systems with minimum of 128 MB RAM and USB port.

Applications

The **e-corder** 821 is an ideal data recorder for use with multiple signal sources. Absolutely no programming is required and the powerful on-line and off-line analysis functions allow results to be quickly extracted from the recorded data. The system can be used to replace paper based chart recorders and data acquisition cards in applications such as Materials Testing, Engineering, Vibration Analysis, Sensor Research and Electrochemistry.


Software


ChartTM and ScopeTM software is shipped with every **e-corder** for both Windows and Macintosh.


Chart TM – lets you use you computer as a multichannel strip chart recorder, polygraph and digital voltmeter.

Scope TM – lets you use your computer as a two-channel storage oscilloscope, or XY plotter.

Specialist software and analysis extensions are also available for use with **e-corder** units.

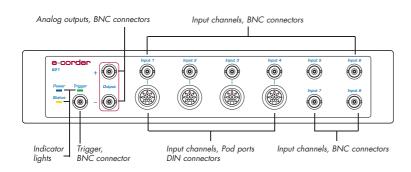

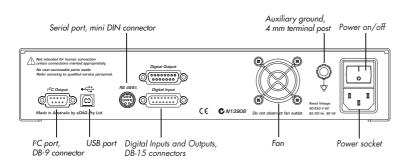

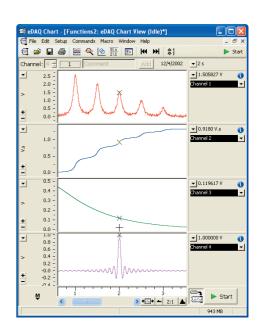

Main Chart window (Windows XP). All software has a user friendly interface and does not require additional programming to develop applications.

Specifications

Number of Input channels:	8
Input configuration:	4 x single-ended BNC, or differential DIN 4 x single-ended BNC
Amplification range:	±2 mV to ±10 V full scale in 12 steps
	Range Resolution (μV) Gain ±10 V 312.5 1 ±5 V 156.25 2 ±2 V 62.5 5 ±1 V 31.25 10 ±0.5 V 15.625 20 ±0.2 V 6.25 50 ±0.1 V 3.125 100 ±50 mV 1.56 200 ±20 mV 0.625 500 ±10 mV 0.3125 1000 ±5 mV 0.15625 2000 ±2 mV 0.0625 5000
Maximum input voltage:	±15 V
Input impedance:	~1 MΩ 47 pF @ DC
Low-pass filtering:	20 kHz fixed 2nd order
AC coupling:	(further filtering via software) DC or 0.15 Hz (software-selectable)
Frequency response (–3 dB):	25 kHz (100 mV - 10 V ranges) 10 kHz (5 mV - 50 mV ranges)
DC drift:	Software corrected zero
CMRR (differential):	>96 dB @100 Hz on 10 mV range
Channel crosstalk:	> -82 dB (typical) between adjacent channel
Input noise:	$<2.4~\mu V_{rms} + 250~\mu _{Vrms}$ /Gain RTI
Pod Connectors	
Connector type:	8-pin DIN
Supply voltage:	±5 V regulated
Maximum current:	50 mA per Pod port
Communications:	2 wire I ² C
Signal input:	Differential or single-ended analog inputs
Sampling	
ADC resolution:	16 bit
Linearity error:	±2 LSB (from 0 to 70 °C)
Maximum sampling rates:	200 kHz for 1 - 2 channels 100 kHz for 3 - 4 channels 40 kHz for 5 - 8 channels
Available sampling rates:	0.003 Hz to 200 kHz using Chart for Mac 2 Hz to 200 kHz using Scope 1 Hz to 200 kHz using Chart for Windows
Output Amplifier	
Output configuration:	Bipolar or dual single-ended
Output resolution:	16 bits

Maximum output current:	2 x ±100 mA continuous
Output impedance:	0.1 Ω typical
Slew rate:	10 V/µs
Settling time:	2 μs (to 0.1% of FSR)
Output range:	±200 mV to ±10 V (software selectable) Range (V) Resolution (μV) ±10 312.5 ±5 156.5 ±2 62.5 ±1 31.25 ±0.5 15.625 ±0.2 6.25
Output Zero error	<0.5 mV
External Trigger	
Trigger mode:	TTL level or contact closure, software selectable
Trigger threshold:	+1.3 V, 50 kΩ input impedance
Hysteresis:	0.3 V
Maximum input voltage:	±12 V
Minimum detectable event:	5 μs
Microprocessor and Data Con	nmunication
CPU:	PowerPC 405 GPr (240 MHz)
RAM:	16 MB DRAM
ROM:	4 Mbit FLASH ROM
Data communication:	USB 2.0 or 1.1 compliant (up to 480 Mbits/s
Expansion Ports	
I ² C expansion port:	Power and control bus for eDAQ Amps (maximum of 500 mA). Interface communications rate of up to 10 000 bits/s.
RS485 Serial Port	RS485 port for factory diagnostic purposes only.
Digital input	8 independent, TTL input lines. 1.2 V threshold, 10 k Ω input impedance. 5 V maximum.
Digital output	8 independent, TTL output level. 8 mA maximum load per line.
Physical Configuration	
Dimensions (w x h x d):	300 x 60 x 300 mm (11.8" x 2.4" x 11.8")
Weight:	4.8 kg (10 lb 8 oz)
Operating voltage:	90 to 250 V AC 50/60 Hz
Nominal power needs:	17 – 30 VA (dependant on number of eDAQ Amps and Pods attached).
	0 to 40 °C

 $\label{eq:WARRANTY: eDAQ Hardware units are supported by a one year warranty} WARRANTY: eDAQ Hardware units are supported by a one year warranty$

